

- For convenient use of up to four Alpha Spectrometers with a single pump
- Screw-type vacuum valves for fine adjustment of pumping speed

The ORTEC Model 576-VM Vacuum Manifold and Control has four valves and terminations for connecting up to four alpha spectrometers to a single pump, in such a way that any one unit may be vented without venting the remainder. The Model 576-VM Vacuum manifold also includes a thermocouple gauge to monitor vacuum. Designed to be installed in a 19-inch rack, it requires 115 V AC, 60 Hz power or 230 V AC, 50 or 60 Hz power.

Specifications

PERFORMANCE

Helium leak checked to $<3 \times 10^{-8}$ atm cm^3/s .

Operating Range: 1 to 5000 millitorr.

CONTROLS

Four screw-type vacuum valves.

INPUTS

Four flexible lines with 1/4 in. Swagelok connectors.

OUTPUT

1-1/8 in. copper tube with NW-25 ISOKF flange termination for connection to user-supplied vacuum pump (minimum pumping speed of 5.6 scfm). May be connected directly to Portable Pump Station. Flange may be cut off with tubing cutter if connection by vacuum rubber hose is desired.

ELECTRICAL AND MECHANICAL



Power Required: 115 V AC, 60 Hz or 230 V AC, 50 or 60 Hz.

Weight

Net: 3 kg (7 lb)

Shipping: 9 kg (20 lb)

Dimensions: 18 cm H x 48 cm W x 53 cm D
(7 x 19 x 21 in.)

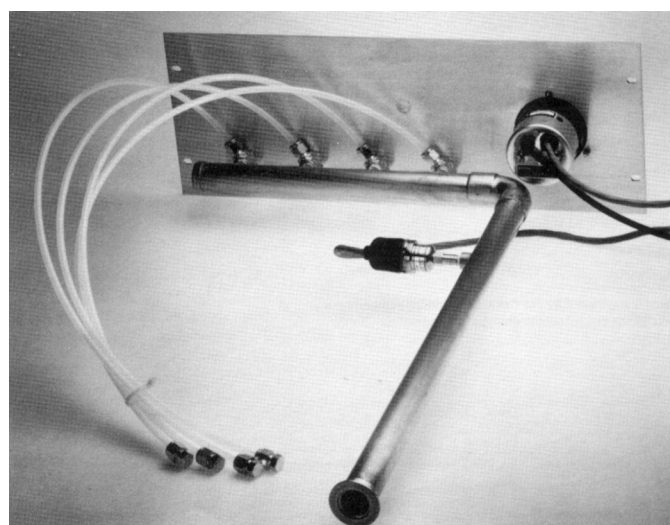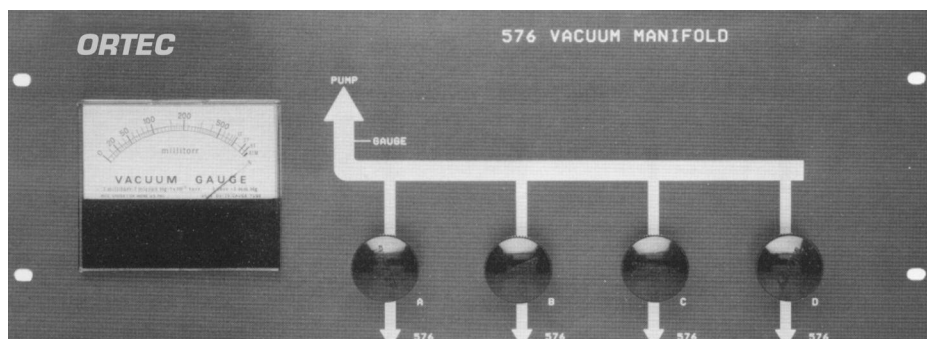

Ordering Information

To order, specify:

Model	Description
576-VM	Vacuum Manifold and Control for 115/230 V AC

Specifications subject to change
100120